

*iDR-L*TM Innovative Direct Radiography



IMAGE

iDR-L | PROVEN DR READER

iDR-L is a scalable design that allows for the highest efficiency **Crystal IP** phosphor screens, 16 Bit available gray scale, and upgradable features that provide users with up to 5 lp/mm maximum resolution.

CAPTURE


XC | ACQUISITION

XC touchscreen acquisition with ICE-4 Enhancement Processing provides all-new features including, "Image Display State" to ensure balanced presentation of both soft tissue, overlapping bone structures, and automatic analysis of image characteristics to optimize processing.

REVIEW

CLARITY PACS

Our fully web-enabled and integrated PACS solutions help transition your practice into a safe, secure, and filmless environment. Clarity PACSTM supports all your current and future imaging needs.

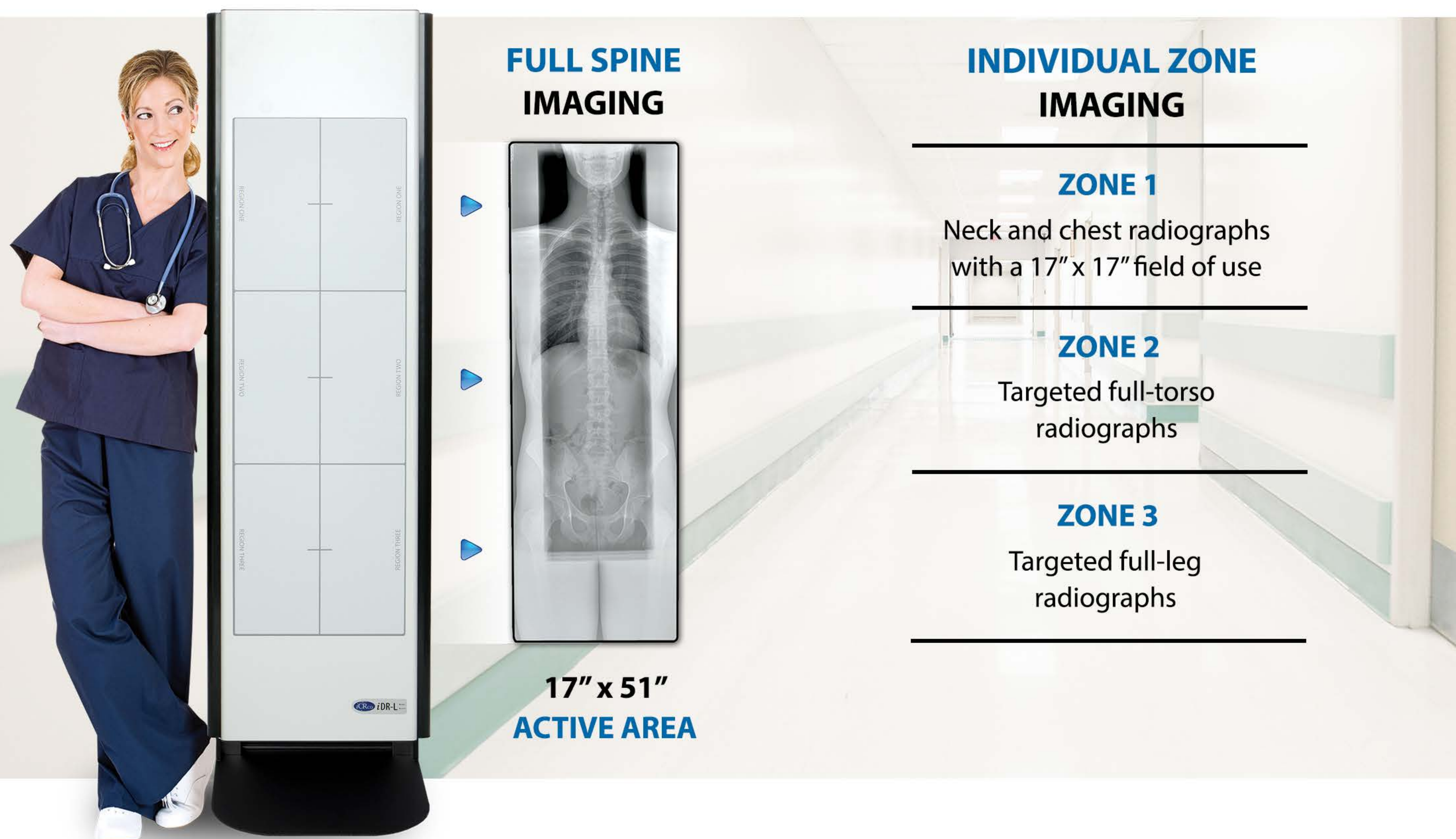

100 Microns

1 Full Spine Detector

3 Regional Scan Areas

16 Bit

Full-Spine Imaging | Higher Resolution | Lower Dose

The exclusive iDR-LTM has a 17"x51" field of use specifically for regional scanning and full-spine imaging studies. Because it takes a single image, you can leverage compensation filters to balance out the exposure between thicker and thinner areas of anatomy.

Physicians can achieve amazing workflow at a much lower cost with its integrated pairing with the robust and feature-rich image processing XCTM software. XC provides superior quality images every time with an option to further manipulate the images.



FULL-SPINE IMAGING IN ONE 17" X 51" EXPOSURE

iDR-L Method: The long-bone radiography platform utilizes a 17" by 51" detector plate that eliminates the need for a cassette. Integrating phosphor plate technology, the iDR-L boasts a higher resolution than standard DR while maintaining the same workflow. This solution sends complete image information directly to the scan processor skipping lengthy image stitching processes. The optional wall stand includes an "EZ-Glide" motorized lift that allows the iDR-L to ascend or descend to any desired position to provide imaging flexibility for chest, abdomen, whole spine, and simultaneous poster anterior (PA) and lateral imaging.

Micron Spot Size	Pixels Per mm	Dots Per Inch (DPI)	Line Pairs Per mm
200	5	127	2.5
155	7.5	163	3.1
100	10	254	5

THE IDR-L CAPTURES AT 15 MEGAPIXELS

Grayscale resolution	16 bits/pixel source file, 65536 shades of gray
Scanning	17" x 51" active area. 14" x 51" option
Wired Data I/F	GigE via optional power & communication tether
Configuration	Wall mounted, upright full-spine imaging
Dimensions	iDR 17: W22.5" X D6.32" x H72.25" iDR 14: W19.5" X D6.32" X H72.25"
Weight	Approx. 100 lbs
Power Source	100-240V AC / 2.5A max; 47-63Hz (universal power supply) 240V option available
iSupply	24V portable DC power supply available, 270 scans on single charge
Heat Generation	Standby 230W, Maximum 1610W
Operating Conditions	Temperature; 0-40°C/32-105°F, temperature change: 0.5°C/min, humidity: 15%-95% RH, magnetic fields: max 1260 µT (in conformance with EN 61000-4-8: Level 3), 10 A/m



Components are made from 98% recyclable parts • iCRco is an ISO 13485 certified company
U.S. and international patents granted • Additional patents pending FDA 510 (K) Cleared
Medical CE mark CE 0086

© 2020 iCRco. All rights reserved. "iDR-L" and "XC" are registered trademarks of iCRco | BR062415 | Design & Specifications are subject to change without notice. Processing and display time dependent on processor speed, RAM disk access time, and video card.

XC™ Acquisition Software Features

- XC - Intuitive touchscreen acquisition
- ICE- 4 Processing - Automated image characteristics analysis for maximum image enhancement
- Image Display State: Automated enhancement of image display at the point of acquisition
- Smart search, sort, and filter options
- Integration with front office management systems like RIS and EMR
- Full set of annotation and measurement tools
- User-preferred settings and privileges

Outer Dimensions



Also available with a
Motorized Wall Stand

Detector available in
14"x51" or 17"x51"